

## **FORMULATION AND STABILIZATION OF ASPIRIN MINI-TABLETS WITH THE AID OF WEAK ACID AND MOISTURE PROTECTIVE COATING**

**\*Raja Subburayalu<sup>1</sup>, Dr. Janakiraman Kunchithapatham<sup>1</sup>, Dr. Ramkumar Pillappan<sup>2</sup>,  
Devi Raja<sup>3</sup>**

<sup>1</sup>Department of Pharmacy, Annamalai University, Annamalai Nagar, Chidambaram, Tamil Nadu, India,

<sup>2</sup>Orchid Chemicals & Pharmaceuticals Ltd., SIPCOT Industrial Park, Irungattukottai, Sriperumbudur 602105, Tamil Nadu, India,

<sup>3</sup>Jaya college of Pharmacy, Thirunindravur, Chennai, Tamil Nadu, India.

Article Received on  
22 August 2013,

Revised on 26 Sept. 2013,  
Accepted on 30 October 2013

**\*Correspondence for  
Author:**

**Raja Subburayalu**  
Department of Pharmacy,  
Annamalai University,  
Annamalai Nagar,  
Chidambaram, Tamil Nadu,  
India  
[rajdevmp@gmail.com](mailto:rajdevmp@gmail.com)

### **ABSTRACT**

The main objective of the present study is to formulate a stable aspirin mini tablets with the aid of a weak acid in core tablet & coated with non aqueous moisture protective coating. Aspirin is highly unstable at alkaline environment & and prone to undergo degradation by hydrolysis. The preformulation study reveals, Aspirin is incompatible with alkali salts, and aspirin itself degrades at exposure condition of elevated temperature and humidity. Whereas, with weak acids, the stability of aspirin is comparatively better. Hence, a weak acid is selected in core tablet. To control the impact of humidity on degradation, a moisture protective layer is coated on core tablet with non aqueous solvent, using conventional coating pan. The coated mini tablets are encapsulated in a hard gelatin capsule shell. The filled

capsules are evaluated for description, assay, dissolution, water by KF at initial and 3months accelerated condition ( $40 \pm 2^\circ\text{C}/75 \pm 5\%\text{RH}$ ), and to conclude the quantity of weak acid & percentage of moisture protective coating required to stabilize the formulation. The formulation with 4mg & 6mg/unit of alginic acid/unit and 4% film coating are failed in stability. The formulation with 8mg/unit of alginic acid & 6% moisture protective coating was found to be stable, and the degradation was controlled.

**Keywords:** Dissolution, Aspirin, mini-tablets, alginic acid, moisture protective coating.

## INTRODUCTION

Aspirin (Acetylsalicylic acid) belong to class of drug referred to as Non-steroidal anti-inflammatory drugs, it's clinically useful as analgesics, anti-inflammation, antipyretic, anti-thrombolytic and anti-rheumatic<sup>[1]</sup>. Aspirin also has an antiplatelet effect by inhibiting the production of thromboxane, which under normal circumstances binds platelet molecules together to create a patch over damaged walls of blood vessels. Because the platelet patch can become too large and also block blood flow, locally and downstream, aspirin is also used long-term, at low doses, to prevent heart attacks, strokes, and blood clot formation in people at high risk for developing blood clots<sup>[2]</sup>. It has also been established that low doses of aspirin may be given immediately after a heart attack to reduce the risk of another heart attack or of the death of cardiac tissue<sup>[3]</sup>.

Aspirin undergoes degradation by hydrolysis. The hydrolysis is catalyzed mainly by alkali, heat in presence of moisture. Aspirin is more stable at pH 3.7 (pKa value of aspirin = 3.7). Degradation of aspirin is aimed to controlled by maintain the product's pH value in line with pKa value of aspirin, with the aid of weak acids. Pellets and mini tablets have gained importance over the years due to their distinctive advantages in both technological and therapeutic aspects. Aspirin In case of aspirin<sup>[4,5]</sup>, direct compression technique has been employed to compress the tablet, because the powder is highly moisture sensitive.

## MATERIALS AND METHODS

### Materials

The following chemicals were obtained from commercial suppliers and used as received: Aspirin (Rhodine, Thailand), Alginic acid (kelacid)(FMC, Europe), Opadry AMB white 80W68912 (colorcon, India), Silicon dioxide (syloid 244 FP)(Grace division, USA), Microcrystalline cellulose (Avicel PH 112) (FMC, USA), Lactose Anhydrous (Supertab 21AN) (DFE pharma, USA), Stearic acid (Merck, Europe), Talc (Luzenac, Italy), pregeatinised starch (Starch 1500, Dow, USA), Isopropyl alcohol and methylene chloride was procured from RFCL Limited., New Delhi, India. All chemicals were reagent grade or higher.

Digital weighing balance (C-220) (make: Sartorius), Mechanical sifter with the screens of ASTM 40# & ASTM 60#, Octagonal blender 4L (Sams tech, India), 16 station rotary

compression machine (Cadmach, India), conventional coating pan(a Remi mechanical propellant stirrer (RA124) (make:Remi), Manual capsule filling machine (MAC 300) (make: Pam machineries), Tray drier (make : Ganson engg), double beam UV Visible spectrophotometer (make: schimadzu), Dissolution test apparatus (Electrolab),

## METHODS

### a. Drug-Excipient polymer compatibility study

Aspirin is individually mixed with different excipient, sifted through ASTM 40#, loaded in to  $40\pm 2^{\circ}\text{C}/75\pm 5\%\text{RH}$  accelerated chamber and exposed for four weeks. Samples are withdrawn after 2 weeks and 4 weeks. The physical admixture of the samples exposed are evaluated at Initial, 2 weeks and 4 weeks at exposed condition for Description & Assay

**Table-1: Aspirin and excipient Compatibility study**

S.No	Drug + Excipient	Ratio
1	Aspirin (D)	-
2	D + Microcrystalline cellulose (Avicel PH112)	1 : 1
3	D + Lactose anhydrous (Supertab 21AN)	1 : 1
4	D + Corn starch (PURE DENT B700)	1 : 0.5
5	D + Pregelatinised starch (Starch 1500)	1 : 0.5
6	D + Alginic acid (Kelacid)	1 : 0.1
7	D + Silicon dioxide (Syloid 244 FP )	1 : 0.5
8	D + Talc	1 : 0.2
9	D + Stearic acid	1 : 0.1
10	D + Opadry AMB white 80W50612	1 : 0.2

### b. Preparation of Mini-Tablets <sup>[6]</sup>,

Formulation of Aspirin mini tablets involves 3 stages

- Stage – I : Blending
- Stage - II : Compression
- Stage - III : Film coating

#### Stage –I Blending

- Aspirin was co-sifted with Advice and lactose through ASTM 40#.

- Alginic acid was co-sifted with pregel starch through ASTM 60#.
- Syloid and Talc was sifted through ASTM 60#.
- The above sifted materials was loaded in octagonal blender and mixed for 15 minutes at 15 rpm.
- Stearic acid was sifted through ASTM 60# and loaded in octagonal blender and mixed for 5 minutes at 15 rpm.
- Each formula was having the batch size of 4000 units.

### **Stage –II Compression**

- The lubricated blend was compressed using 2.5mm multi-tip punch with the target weight of 20mg/unit and 8 units/unit dose.
- The compressed mini-tablets were evaluated for hardness, thickness and disintegration time.
- The lot size for barrier coating was 4000 units.

### **Stage-III: Film coating**

- Opadry AMB white 80W50612 suspended in Isopropyl alcohol and methylene chloride admixture under stirring.
- Stirring for continued for 30 minutes.
- The resultant suspension was coated on compressed mini-tablets with different percentage weight gain by using conventional coating pan.
- During the preparation of coating suspension the 10% of excess was prepared to recover the loss during practical work. And the coating suspension was compressed mini-tablets using conventional coating pan until weight gain was achieved and % yield was calculated.
- The solid content of film coating suspension was 7% w/w

### **Encapsulation**

The coated mini-tablets were filled in to size “0” hard gelatin capsules, and evaluated for assay and dissolution.

Table-2 Composition of Aspirin Film coated Mini-tablets

S.No	Ingredients	F1	F2	F3	F4
1	Aspirin (Rhodine 2080A)	75	75	75	75
2	Corn starch NF	12	12	12	12
3	Alginic acid	4	6	8	10
4	Lactose anhydrous (Supertab 21 AN)	25	23	21	19
5	Microcrystalline cellulose (Avicel PH112)	38	38	38	38
6	Silicon dioxide (Syloid 244FP)	2	2	2	2
7	Talc	2	2	2	2
8	Stearic acid	2	2	2	2
	<b>Sub total</b>	<b>160</b>	<b>160</b>	<b>160</b>	<b>160</b>
<b>Film coating</b>					
S.No	Ingredients	F-1A	F-2A	F-3A	F-4A
1	Opadry AMB white	6.4	6.4	6.4	6.4
2	Methylene chloride	qs	qs	qs	qs
3	Isopropyl alcohol	qs	qs	qs	qs
	<b>Sub total</b>	<b>166.4</b>	<b>166.4</b>	<b>166.4</b>	<b>166.4</b>
S.No	Ingredients	F-1B	F-2B	F-3B	F-4B
1	Opadry AMB white	9.6	9.6	9.6	9.6
2	Methylene chloride	qs	qs	qs	qs
3	Isopropyl alcohol	qs	qs	qs	qs
	<b>Total</b>	<b>169.6</b>	<b>169.6</b>	<b>169.6</b>	<b>169.6</b>
S.No	Ingredients	F-1C	F-2C	F-3C	F-4C
1	Opadry AMB white	12.8	12.8	12.8	12.8
2	Methylene chloride	qs	qs	qs	qs
3	Isopropyl alcohol	qs	qs	qs	qs
	<b>Total</b>	<b>172.8</b>	<b>172.8</b>	<b>172.8</b>	<b>172.8</b>

Table-3: In-Process parameters at various steps

Parameters	Compression	Coating
<b>Description</b>	White to off-white circular biconvex mini tablets	White to off-white circular biconvex mini tablets
<b>Uniformity of tablet weight</b>	20mg $\pm$ 2mg	Target weight $\pm$ 10%
<b>Hardness (N)</b>	10-20N	20-40N
<b>Thickness (mm)</b>	3.1- 3.5mm	3.1- 3.5mm
<b>Disintegration time</b>	NMT 15 min	NMT 15 min

### Stability study<sup>7</sup>

Stability testing of drug products begins as a part of drug discovery and ends with the emise of the compound or commercial product. FDA and ICH specify the guidelines for stability testing of new drug products, as a technical requirement for the registration of pharmaceuticals for human use. The ICH Guidelines have established different temperatures and period of stability testing. The coated mini-tablets of formulation F1A to F4C were filled in size '1' capsules, packed in HDPE bottle, and loaded on stability chamber as per ICH guidelines, as mentioned in table-4.

**Table-4: ICH guidelines for Stability Study**

Study Storage	Condition	Time
Long term	25°C±2°C / 60% RH±5% RH	12 month
Intermediate	30°C±2°C /65% RH±5% RH	12 months
Accelerated	40°C±2°C /75% RH±5% RH	6 months

To evaluate the impact of weak acid and moisture protective layer, the product is evaluated at accelerated stability condition.

### Evaluation of Aspirin film coated mini tablets

Uniformity of Tablet Weight Test <sup>[8]</sup>: Ten capsules from the batch were randomly selected, individual weight of the selected representative was determined using a digital electronic balance. The average tablet weight and the standard deviation from the mean were calculated.

Capsule disintegration Test <sup>[8]</sup>: Six capsules randomly selected were introduced into the six baskets of the disintegration testing apparatus (Electrolab, India). The disintegrating medium was de ionized water maintained at 37°C + 1.0°C. The time taken for each capsule to disintegrate to break up into a smaller units and passes through the screen mesh orifices at the bottom of the basket was recorded.

### Dissolution <sup>[9, 10]</sup>

Acetate buffer at pH 4.5 was used as the dissolution medium as specified in the British Pharmacopoeia. It was prepared by mixing 29.9 g of sodium acetate with 16.6 mL of glacial acetic acid and sufficient distilled water to produce 10 L. Sodium acetate and glacial acetic acid were analytical grade and purchased from Sigma–Aldrich UK. Pure aspirin powder (Acetylsalicylic Acid BP, Sigma– Aldrich UK) was dissolved in the acetate buffer to make a

series of standard calibration solutions with different concentrations for development of a calibration curve using a UV spectrophotometer at 265 nm.

**Dissolution Testing** In vitro dissolution was carried out via USP Apparatus 2 (paddle) at a speed of 75 rpm in 900 ml of dissolution medium (pH 4.5 acetate buffer) maintained at  $37 \pm 0.5$  °C using a water bath fitted with a variable- speed stirrer and heater (Erweka DT6). Selection of 75- rpm rotation speed was based on the British Pharmacopoeia guideline (22). Samples ( $5 \pm 0.1$  mL) were taken manually at 10, 20, 30, 45, 60, and 90 min and replaced with an equal volume of fresh medium to maintain a constant dissolution volume. The samples were filtered, and the absorbance was measured at 265 nm using a UV spectrophotometer (PU 8625 UV/VIS spectrophotometer). The drug concentration determined by the calibration model was used to calculate the total mass of the drug released in the medium. In this work, the dissolution profiles are represented as the cumulative percentages of the amount of drug released at each sampling interval. Each profile is the average of six individual tablets.

### Assay <sup>[8]</sup>

Aspirin film coated mini-tablets from the capsule were dispersed in to 190 ml of pH 6.8 phosphate buffers by ultrasonication for 30 minutes followed by 10 minutes stirring using magnetic stirrer. The solution was then filtered and the residues over filter paper were washed with 10 ml phosphate buffer. The solution was then diluted up to suitable concentration and absorbance was measured using double beam UV-VIS Spectrophotometer at 289 nm.

## RESULTS & DISCUSSIONS

### Drug-Excipient compatibility study

The Drug – Excipient compatibility study is conducted for assay and description, the results are tabulated in **Table 5**.

Drug : Excipient	Ratio	Initial		2 Week 40 $\pm$ 2°C/75 $\pm$ 5%RH		4 Week 40/75 40 $\pm$ 2°C/75 $\pm$ 5%RH	
		Description	Assay	Description	Assay	Description	Assay
Aspirin (D)	-	White to off white powder	99.6	Off white powder	98.7	Off white to pale pink powder	93.5

D + Microcrystalline cellulose (Avicel PH112)	1 : 1	Off - white powder	101.1	pale pink colored powder	95.8	pale pink colored powder	91.2
D + Lactose anhydrous (Supertab 21AN)	1 : 1	Off - white powder	99.9	pale pink colored powder	96.5	pale pink colored powder	92.1
D + Corn starch (PURE DENT B700)	1 : 0.5	Off - white powder	101.1	pale pink colored powder	95.8	pale pink colored powder	92.1
D + Pregealtinised starch (Starch 1500)	1 : 0.5	Off - white powder	99.6	pale pink colored powder	98.5	pale pink colored powder	93.6
D + Alginic acid (Kelacid)	1 : 0.2	Off - white powder	100.2	Off white powder	99.8	Off white powder	99.7
D + Silicon dioxide (Syloid 244 FP )	1 : 0.5	Off - white powder	99.6	pale pink colored powder	97.5	pale pink colored powder	94.3
D + Talc	1 : 0.2	Off - white powder	99.9	Off white powder	97.6	Off white to pale pink powder	95.2
D + Stearic acid	1 : 0.1	Off - white powder	100.2	Off white powder	97.5	Off white to pale pink powder	95.3
D + Opadry AMB white	1: 0.2	Off - white powder	100.3	Off white powder	96.5	Off white to pale pink powder	94.3

The results of compatibility study reveal that the drug is not stable at exposed condition as alone and with excipient. But the same was stabilized using alginic acid. Around 5-8% drop in potency is observed after 4 week direct exposure. Hence, a moisture protective layer and alginic acid is required to stabilize aspirin.



## Physical Characterization of Aspirin mini-tablets filled in capsules

Table 6: Physical Characterization of Aspirin mini-tablets filled in capsules

Parameters	Specification	Batch Number			
		F1A	F2A	F3A	F4A
Number of mini-tablets per capsule	8	8	8	8	8
Uniformity of tablet weight	Mean $\pm$ SD	20.9 $\pm$ 0.10	20.8 $\pm$ 0.08	20.8 $\pm$ 0.12	20.8 $\pm$ 0.09
Disintegration time	NMT 15 min	4min 20 sec	3 min 45 sec	3 min 35 sec	3 min 15 sec
Parameters	Specification	Batch Number			
		F1B	F2B	F3B	F4B
Number of mini-tablets per capsule	8	8	8	8	8
Uniformity of tablet weight	Mean $\pm$ SD	21.2 $\pm$ 0.13	21.3 $\pm$ 0.08	21.2 $\pm$ 0.10	21.2 $\pm$ 0.07
Disintegration time	NMT 15 min	5min 5sec	4 min 35 sec	4min 25sec	4min 45 sec
Parameters	Specification	Batch Number			
		F1C	F2C	F3C	F4C
Number of mini-tablets per capsule	8	8	8	8	8
Uniformity of tablet weight	Mean $\pm$ SD	21.7 $\pm$ 0.10	21.7 $\pm$ 0.07	21.6 $\pm$ 0.11	21.7 $\pm$ 0.07
Disintegration time	NMT 15 min	5min 30 sec	4min 55 sec	4min 50 sec	5min 10sec

**Physical Description of mini-tablets filled in capsules****Table 7: Comparison of physical description of Aspirin Mini-tablets filled in capsules 75mg (Initial Vs 3 months  $40 \pm 2^\circ\text{C}/75 \pm 5\% \text{ RH}$ )**

Batch Number↓	Physical Description	
	Initial	3 months $40 \pm 2^\circ\text{C}/75 \pm 5\% \text{ RH}$
ASCL-C-F1A (4mg alginic acid + 4% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	pale pink colored mini tablets filled in size"0" hard gelatin capsules
ASCL-C-F1B (4mg alginic acid + 6% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	Pale pink colored mini tablets filled in size"0" hard gelatin capsules
ASCL-C-F1C (4mg alginic acid + 8% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	pale colored mini tablets filled in size"0" hard gelatin capsules

Batch Number↓	Physical Description	
	Initial	3 months $40 \pm 2^\circ\text{C}/75 \pm 5\% \text{ RH}$
ASCL-C-F2A (6mg alginic acid + 4% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	pale colored mini tablets filled in size"0" hard gelatin capsules
ASCL-C-F2B (6mg alginic acid + 6% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	Pale colored mini tablets filled in size"0" hard gelatin capsules
ASCL-C-F2C (6mg alginic acid + 8% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	Pale colored mini tablets filled in size"0" hard gelatin capsules

Batch Number↓	Physical Description	
	Initial	3 months $40 \pm 2^\circ\text{C}/75 \pm 5\%$ RH
ASCL-C-F3A (8mg alginic acid + 4% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	White to off white colored mini tablets filled in size"0" hard gelatin capsules
ASCL-C-F3B (8mg alginic acid + 6% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	White to off white colored mini tablets filled in size"0" hard gelatin capsules
ASCL-C-F3C (8mg alginic acid + 8% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	White to off white colored mini tablets filled in size"0" hard gelatin capsules

Batch Number↓	Physical Description	
	Initial	3 months $40 \pm 2^\circ\text{C}/75 \pm 5\%$ RH
ASCL-C-F4A (10mg alginic acid + 4% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	White to off white colored mini tablets filled in size"0" hard gelatin capsules
ASCL-C-F4B (10mg alginic acid + 6% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	White to off white colored mini tablets filled in size"0" hard gelatin capsules
ASCL-C-F4C (10mg alginic acid + 8% film coating)	White to off white colored mini-tablets pellets filled in size"0" hard gelatin capsules	White to off white colored mini tablets filled in size"0" hard gelatin capsules

Results indicated that Batches F3A, F3B, F3C, F4A, F4B & F4C are not having any change in physical description in 3 months accelerated condition. F1A, F1B, F1C, F2A, F2B & F2C were failed in physical description at 3 months accelerated condition. From the results it was concluded that minimum 8mg alginic acid is required to have stable formulation.


### Drug content

Initial and Accelerated stability samples of the filled capsules were evaluated for Drug content & acid resistance. The results are tabulated in Table:

Table 8: Assay, Water by KF and Dissolution: (Initial Vs 3M Accelerated condition)

Batch Number↓	Initial			3M Accelerated condition		
	Assay	Water by KF	Dissolution	Assay	Water by KF	Dissolution
Limit	90-110% of Label claim	NMT 4% w/w	NLT 80%(Q) in 30 min	90-110% of Label claim	NMT 4% w/w	NLT 80%(Q) in 30 min
F1A	100.2 ± 0.36	1.25	93 ± 2.1	90.7 ± 0.67	5.6	62 ± 0.6
F2A	100.4 ± 0.2	1.42	97 ± 2.1	91.5 ± 0.66	4.8	62 ± 0.6
F3A	99.8 ± 0.25	1.5	96 ± 1.5	97.1 ± 0.53	2.8	90 ± 1.2
F4A	99.9 ± 0.35	1.52	96 ± 0.6	98.3 ± 0.23	3.2	85 ± 2.1
F1B	100.7 ± 0.81	1.32	93 ± 3.8	90.2 ± 0.81	3.1	63 ± 1.7
F2B	100.5 ± 0.25	1.25	96 ± 1.5	91.7 ± 0.35	3	76 ± 2.0
F3B	100.5 ± 0.06	1.24	97 ± 1.2	99.2 ± 0.15	1.9	92 ± 1.5
F4B	100.7 ± 0.12	1.25	99 ± 1.0	96.5 ± 0.1	1.85	89 ± 1.0
F1C	100.7 ± 0.5	1.25	96 ± 2.0	92.6 ± 0.12	2.8	69 ± 1.5
F2C	100.3 ± 1.19	1.32	98 ± 1.5	96.7 ± 0.17	2.4	77 ± 1.2
F3C	100.3 ± 0.42	1.22	98 ± 1.5	99.1 ± 0.15	1.85	90 ± 1.2
F4C	100.5 ± 0.06	1.2	96 ± 2.1	98.3 ± 0.2	1.9	87 ± 1.2

\*Listed value indicates mean value of results and Standard deviation (Where n=3)

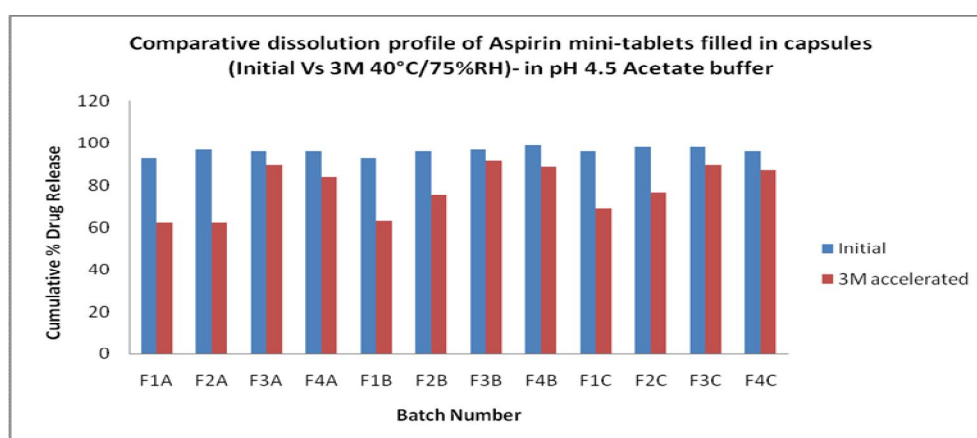

Graph-1: Comparative % drug release of Aspirin mini-tablets filled in capsules Initial Vs 3M 40/75 - in pH 4.5 Acetate buffer, for 30 minutes.

The assay, water content result indicates that batch number F1A, F2A, F1B and F2B are showing high degradation after 3M accelerated condition. Hence 4mg of alginic acid is not efficient to stabilize the formulation of aspirin mini-tablets. Formulation F1A, F2A, F1B, F2B, F1C and F2C are failing in dissolution after 3M accelerated condition. Hence 4 mg and 6 mg of alginic acid are not efficient to stabilize the formulation of aspirin mini-tablets. The formulations of F3A, F3B, F3C (with 8mg of alginic acid ) and F4A, F4B and F4C (10mg of alginic acid) are showing assay, water by KF and dissolution within the specified limit after 3M accelerated condition. The formulation F3B is comparatively better and the drop in dissolution is not observed.

## CONCLUSION

The preformulation result concludes that aspirin is highly unstable at alkaline condition and humidity, and degrades by hydrolysis. Aspirin is highly stable at the pH value of 3.7(pKa value of aspirin). A weak acid was evaluated from 4mg/unit to 10mg/unit. With the aid of weak acid 8mg per unit the product's pH was attained to pH 3.7, and the formulation was found to be stable at accelerated condition at 3month, and significant difference was not observed in comparison to initial. The formulation was evaluated with the moisture protective layer coating of 4% w/w to 8%w/w buildup. Whereas, the formulation with 6%w/w & 8%w/w moisture protective layer is found to be stable. Hence, the formulation with 8mg/unit of alginic acid and a moisture protective layer coating of 6% w/w was finalized as stable formulation.

## REFERENCES

1. Grotta JC, Current Medical and Surgical Therapy for Cerebro-vascular disease, Eng. J. Med., 317, 1987, 1505-16.
2. Lewis HD, Davis JW, Archibald DG, Steinke WE, Smitherman TC, Doherty JE, Schnaper HW, LeWinter MM, Linares E, Pouget JM, Sabharwal SC, Chesler E, DeMots H. The New Eng J Med, 1983; 309 (7): 396– 403.
3. Krumholz, Harlan M, Martha J. Radford, Edward F. Ellerbeck, John Hennen, Thomas P. Meehan, Marcia Petrillo, Yun Wang, Timothy F. Kresowik, Stephen F. Jencks. Circulation, 1995; 92 (10): 2841-7.
4. British Pharmacopoeia, Vol.3, 2009.
5. Physician Desk Reference, 61st ed. 2007: 1869.

6. NOORANA TEHSEEN, , VINAY RAO, MOHD ABDUL HADI DESIGN AND CHARACTERIZATION OF TWICE DAILY MINI-TABLETS FORMULATION OF PREGABALIN *International Journal of Pharmaceutical science*, Voume 5, Supplement 1, page 168-175
7. ICH, GUIDELINES Q1C, "Guidance for industry, stability testing of new dosage form" November 1996. <http://www.ich.org/about/organisation-of-ich/coopgroup/asean/topics-underharmonisation/article/stability-study.html>
8. Oyeniyi and Itiola OA, Pharmaceutical evaluation of direct compressible acetyl salicylic acid tablets containing sawdust microcrystalline cellulose , *International journal of biology, pharmacy and applied science*. April, 2012, 1(3): 195-203
9. Comparison of In Vitro Dissolution Tests for Commercially Available Aspirin Tablets Fawzia Khan, Mingzhong Li\*, and Walkiria Schlindwein School of Pharmacy, De Montfort University, the Gateway, Leicester, LE1 9BH, UK
10. United States Pharmacopoeia. The National Formulary, (USP33/NF28) ed, Rockville: Maryland, The United States Pharmacopoeial Inc., 2010.